

物理实验教学中心

Physics Experiment Center



EQUAL THICKNESS INTERFERENCE

Li Bin
NJUPT
2025 Spring

Purposes:

1. Observation and measurement of the equal thickness interference image.
2. The basic regulation and measurement operation of the microscope.
3. Measure the curvature radius of lens using **Newton ring**.
4. Learn to use graphical method and differential method for data analysis.

Principles

o is the touch point, $e = 0$ at o point.

Optical path difference (OPD): δ

$$\delta = 2e + \frac{\lambda}{2} = 2Re - e^2$$


$$r^2 = R^2 - (R - e)^2 \approx 2Re$$

$$\delta = \frac{r^2}{R} + \frac{\lambda}{2} = \begin{cases} 2k\frac{\lambda}{2} \\ (2k+1)\frac{\lambda}{2} \end{cases}$$

Bright
rings

Dark
rings

Newton ring

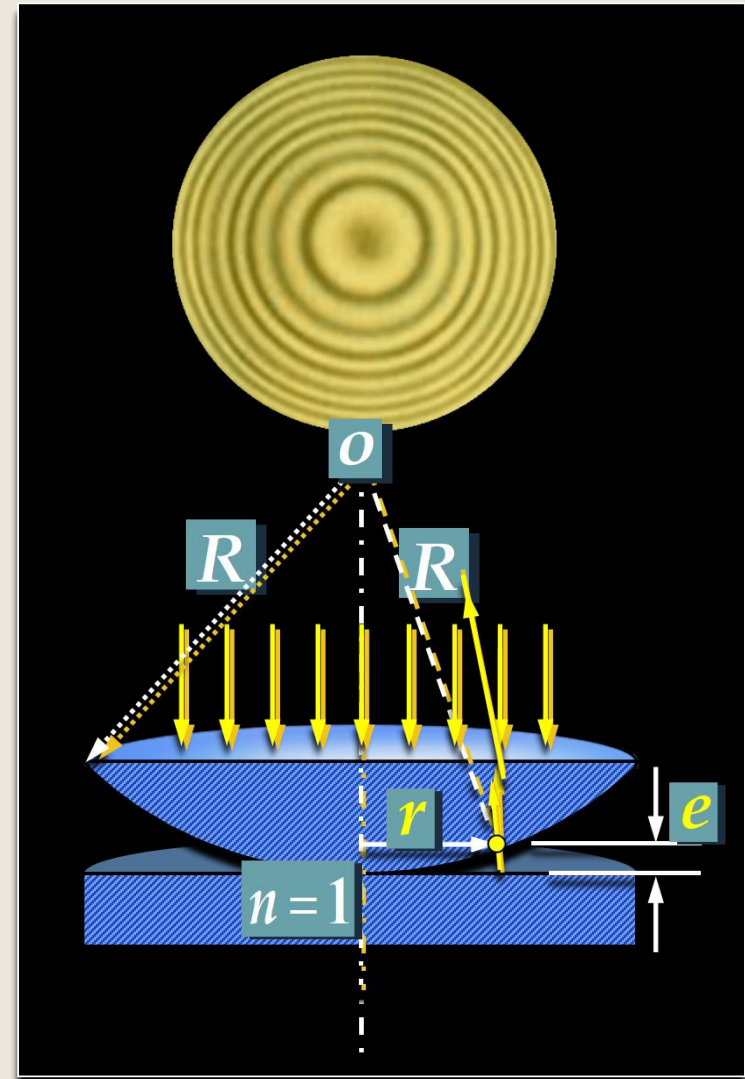

Convex len


Measurement of the curvature radius (R) of convex lens using Newton ring

For Dark rings:

$$D_k^2 = (4R\lambda)k$$

$$R = \frac{D_m^2 - D_n^2}{4(m - n)\lambda}$$



m, n is the order

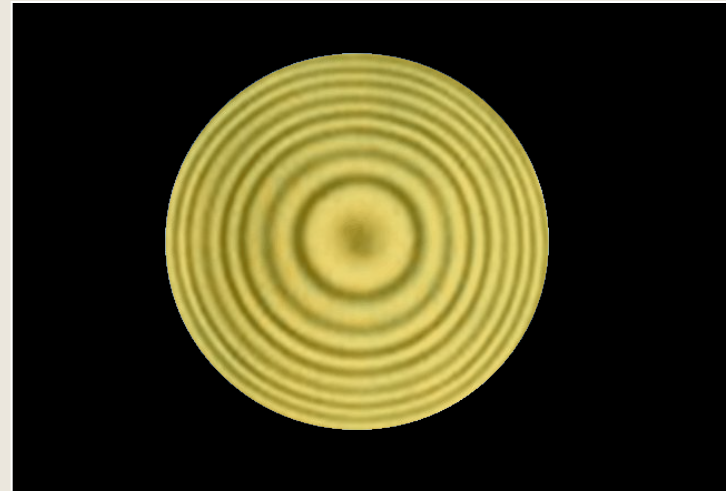


Features: Bigger the level k , more dense the stripes. ($k = 0$ at center)

Instruments

Microscope, Newton Convex lens, Natrium lamp ($\lambda = 589.3\text{nm}$)

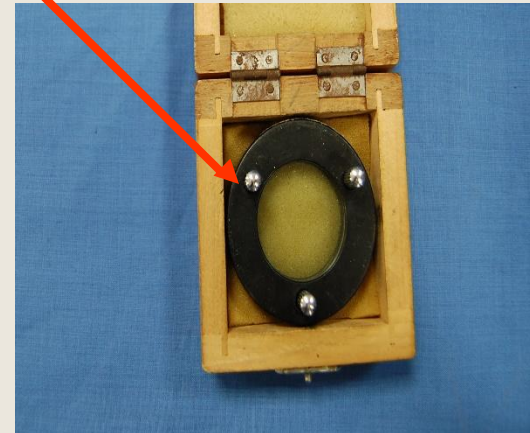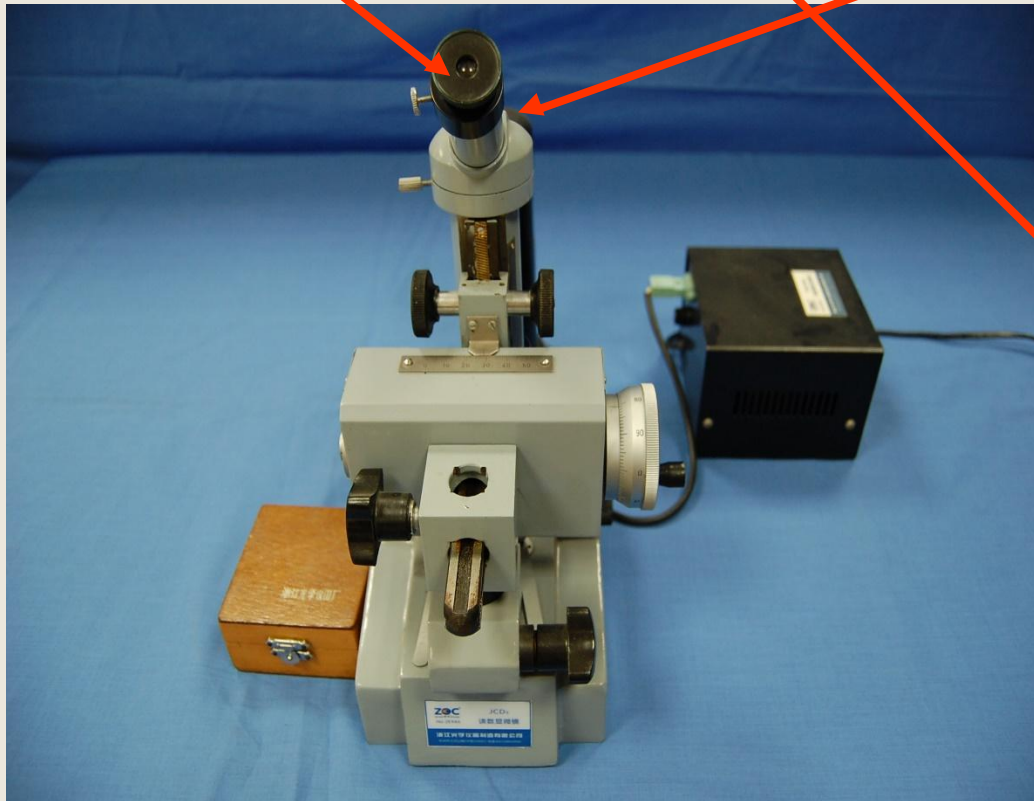

How to use microscope?

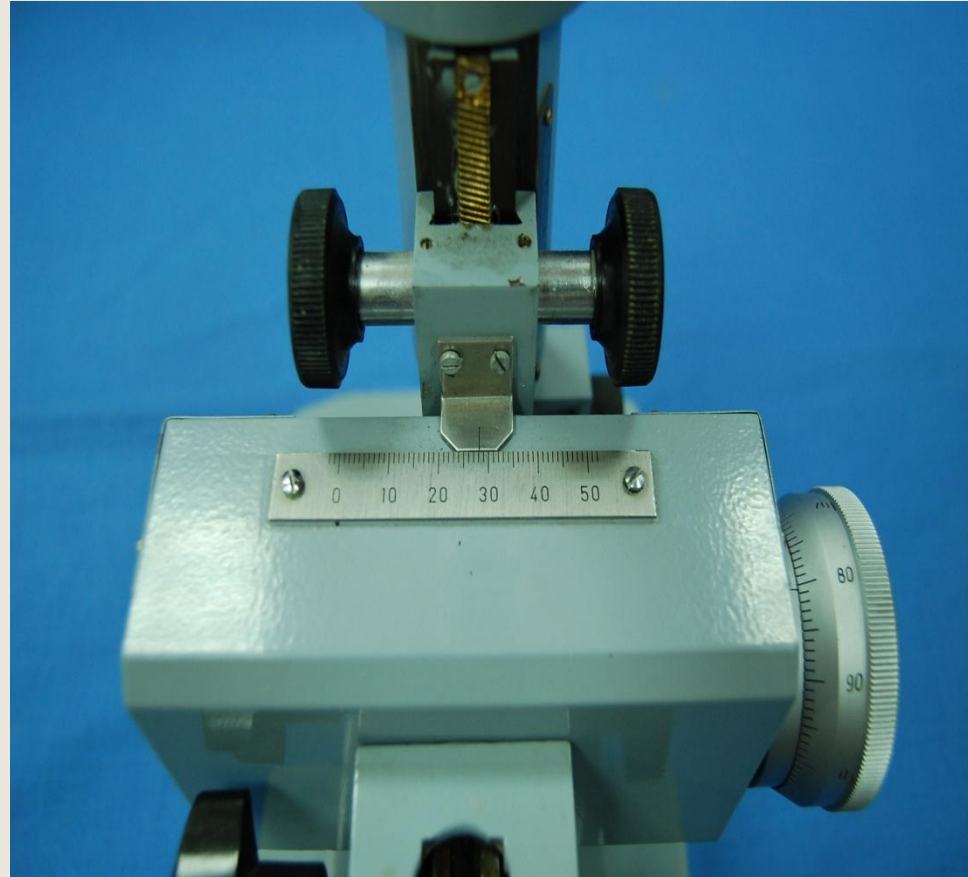





① Corse adjustment

② Fine adjustment

- a. Regulate eyepiece, make the cross wire clear
- b. Rotate transparent mirror, make the view of the microscope most bright.

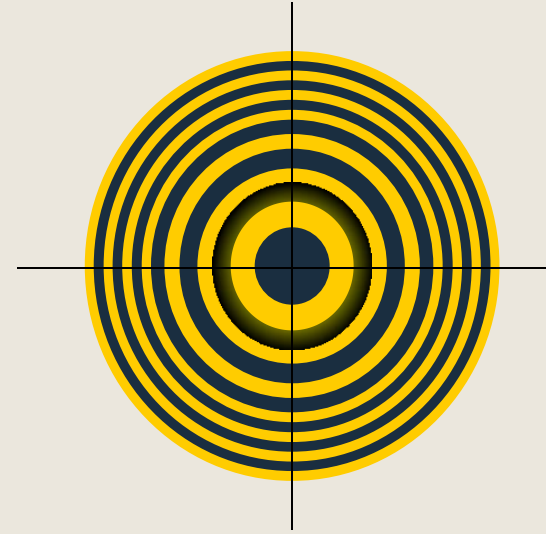

c. Rotating hand wheel from bottom to top slowly until the stripes clear.

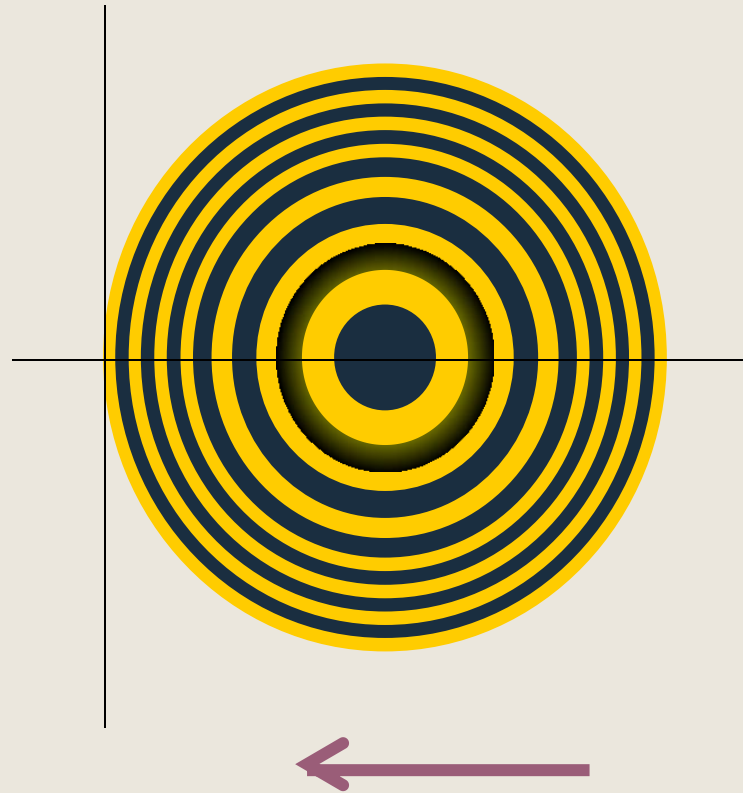

To eliminate the return difference, we firstly move the cross wire to 35th dark ring on the right side, then turn back to the 30th dark ring, start to record the data as X_{30} , X_{25} , X_{20} , X_{15} , X_{10} , X_5 of both sides (Right to Left).

Operation video



Table I

Order K		30	25	20	15	10	5
Data (cm)	Left X_K						
	Right X_K'						
Diameter of rings $D_K= X_K-X_K' $ (cm)		$D_{30}=$	$D_{25}=$	$D_{20}=$	$D_{15}=$	$D_{10}=$	$D_5=$
D_K^2 (cm ²)		$D_{30}^2=$	$D_{25}^2=$	$D_{20}^2=$	$D_{15}^2=$	$D_{10}^2=$	$D_5^2=$
$D_m^2-D_n^2$ (cm ²) ($m-n=10$)		$D_{30}^2-D_{15}^2=$		$D_{25}^2-D_{10}^2=$		$D_{20}^2-D_5^2=$	
$\overline{D_m^2-D_n^2}$ (cm ²)		Average value of $D_m^2-D_n^2$					
\overline{R} (cm)							

$$\overline{R} = \frac{\overline{D_m^2 - D_n^2}}{4(m-n)\lambda} \quad \lambda = 589.3\text{nm}, m-n=15$$

Plotting

$$D_k^2 = (4R\lambda)k$$

Use k as x axis, D_k^2 as y axis, plot $k \sim D_k^2$ graph, get the **slope** of the line: $(4R\lambda)$, then calculate the value of R : $R = \text{slope} / (4\lambda)$. ($\lambda = 589.3\text{nm}$)

Compare this R with the previous one.

Assignment

- Please write a 300-word essay to describe equal thickness interference, complete Table I and plot $k \sim D_k^2$ graph.
- DL: Next class, April 24, 2025

- Useful links:

1. <https://github.com/bliseu/phylab>



END