# 物理实验数学中心

Physics Experiment Center



Website: https://sourceforge.net/projects/phy-njupt/files/2017Spring/

# Oscilloscope

Li Bin


NJUPT


# Experimental Goals


- 1. Adjustment and use of oscilloscope.
- 2. Learn to use oscilloscope to observe voltage waveform.
- 3. Observation of Lissajous figures.

#### GDS-2062 Oscilloscope



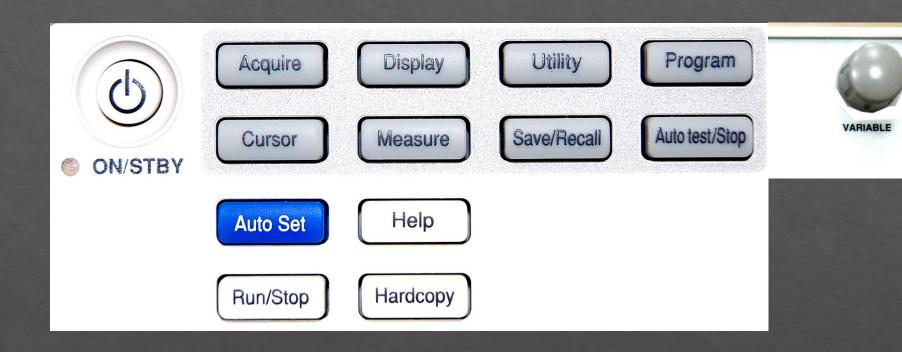

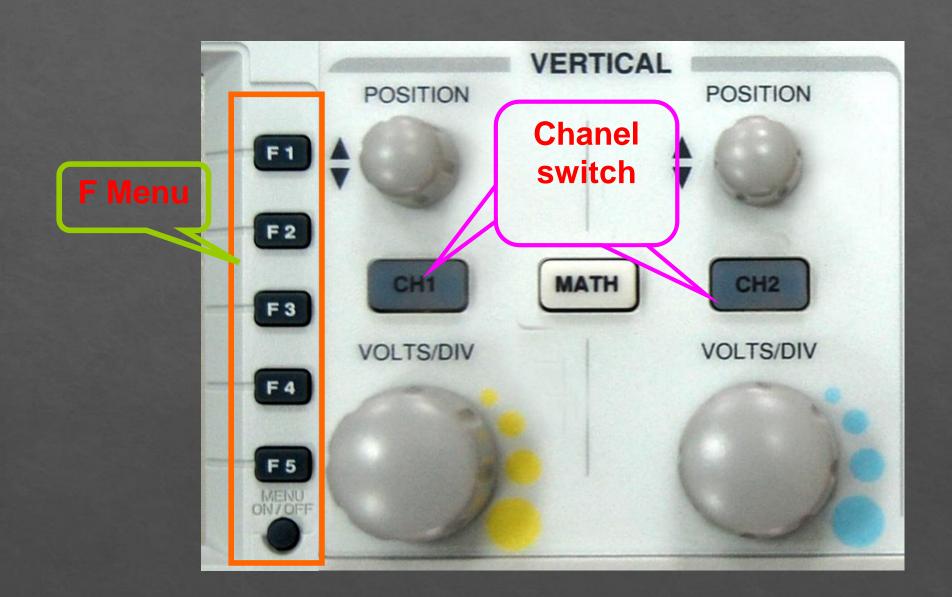


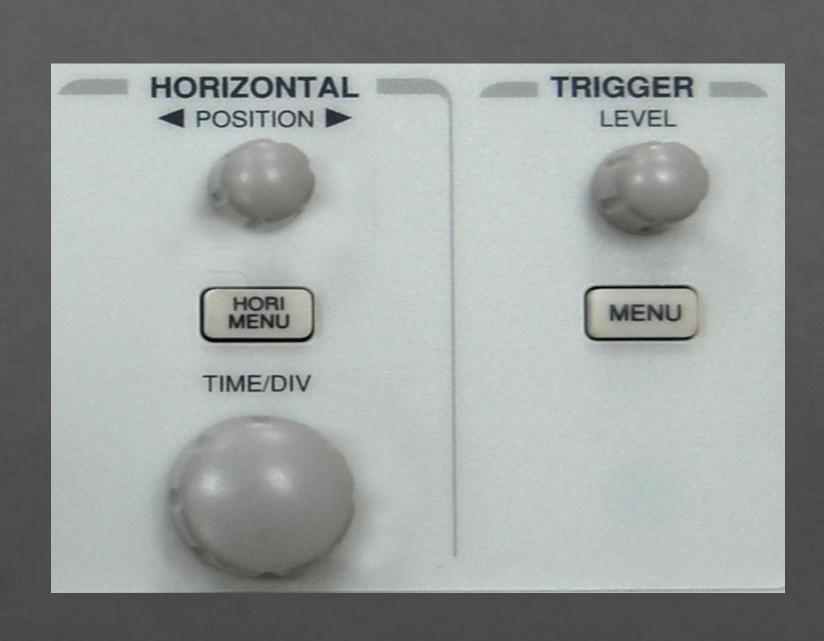

#### F05A Signal generator





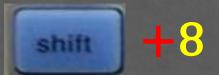

Chanel A:



Chanel B:





# Contents and Steps:

# 1. Settings of Signal generator:


- (1) Set the frequency: + value + unit
- (2) Set the voltage: + value + unit
- (3) Set the phase difference between channel A & B:



# 2. Observe voltage waveform on Oscilloscope

(1)Power on;

(2) Open the channels (corresponding lights on), press Auto Set or adjust the VOLTS/DIV &TIME/DIV buttons, observe the voltage waveform.

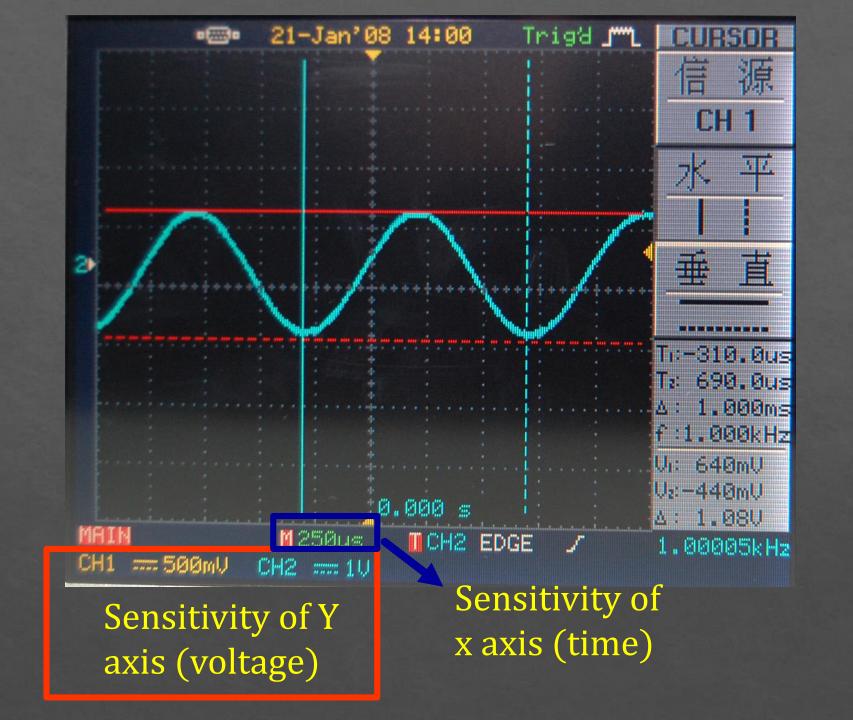

#### 3. Lissajous figures

(1) Turn on two channels;





On signal generator

### Phase differences:


On signal generator



#### Ny=6; Lissajous figures

Nx:Ny=2:3;

Fx:Fy=Ny:Nx=3:2

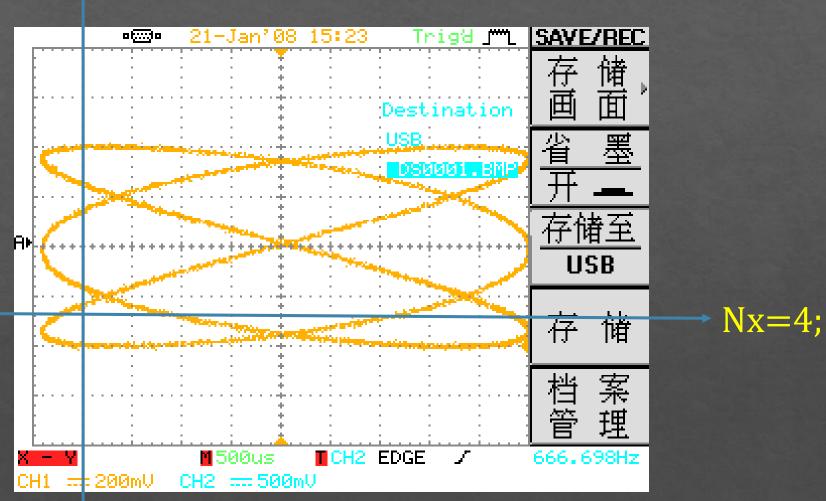

#### TABLES

$$U_p = \frac{1}{2} U_{p-p} = \sqrt{2} U$$

# 1. The voltage of sinusoidal signals

NO.	Voltage (V)	Sensitivity of Y axis (V)	D <sub>y</sub>	U <sub>p-p</sub> (v)	U <sub>p</sub> (v)	U (v)
1	2.5					
2	3.0					
3	4.5					
4	5.0					

## 2. The period of sinusoidal signals

NO.	f (H <sub>Z</sub> )	Sensitivity of X axis (ms)	$D_x$	T (ms)
1	400			
2	2000			
3	8000			
4	15000			

## 3. Lissajous figures

$$f_x = 3 \text{ kHz}, f_y = 6 \text{ kHz}$$

Phase diff.	$0_{\bar{o}}$	90∘	180º	270⁰	360º
Lissajous figures					

### END