

物理实验教学中心

Physics Experiment Center



Dielectric constant measurement

Li Bin

NJUPT


Purposes:

- 1、 Handle the principles and method of the measurement of dielectric constants.
- 2、 Learn the method of data processing.

Principles:


1、 Measure the dielectric constant of solid dielectrics by alternating current bridge

In the air

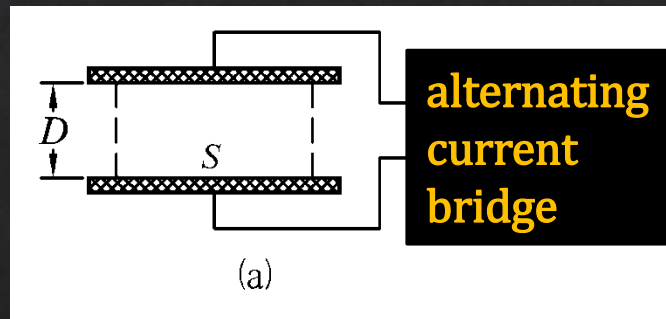

$$C_1 = C_0 + C_{B1} + C_{F1}$$

Insert dielectric sample

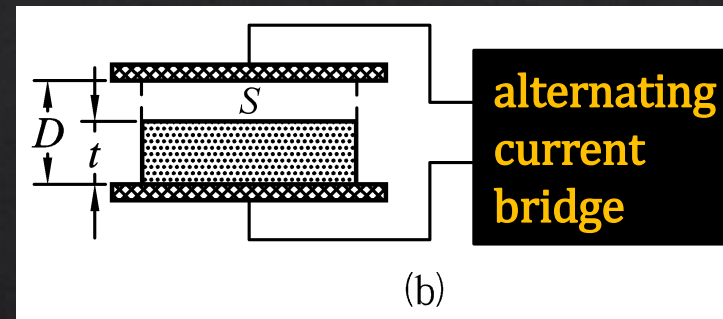

$$C_2 = C_C + C_{B2} + C_{F2}$$

We have these functions:

$$C_0 = \frac{\varepsilon_0 S}{D}$$

$$C_C = \frac{\varepsilon_r \varepsilon_0 S}{t + \varepsilon_r (D-t)}$$

$$C_{B1} = C_{B2}$$

$$C_{F1} = C_{F2}$$

$$C_C = C_2 - C_1 + C_0$$

$$\varepsilon_r = \frac{C_C \cdot t}{\varepsilon_0 S - C_C (D-t)}$$

2、 Calculations of dielectric constant and capacitance:

Vacuum dielectric constant ϵ_0 , S_0 is plate area, D is distance between plates, the capacitance C is shown as:

$$C = \frac{\epsilon_0 S_0}{D} + C_F$$

Set $C = y$, $\frac{1}{D} = x$


$$y = a + bx$$

We got $a = C_F$, $b = \epsilon_0 S_0$

Finally, we have $\epsilon_0 = b/S_0$

Instruments:

1、 Solid electrode

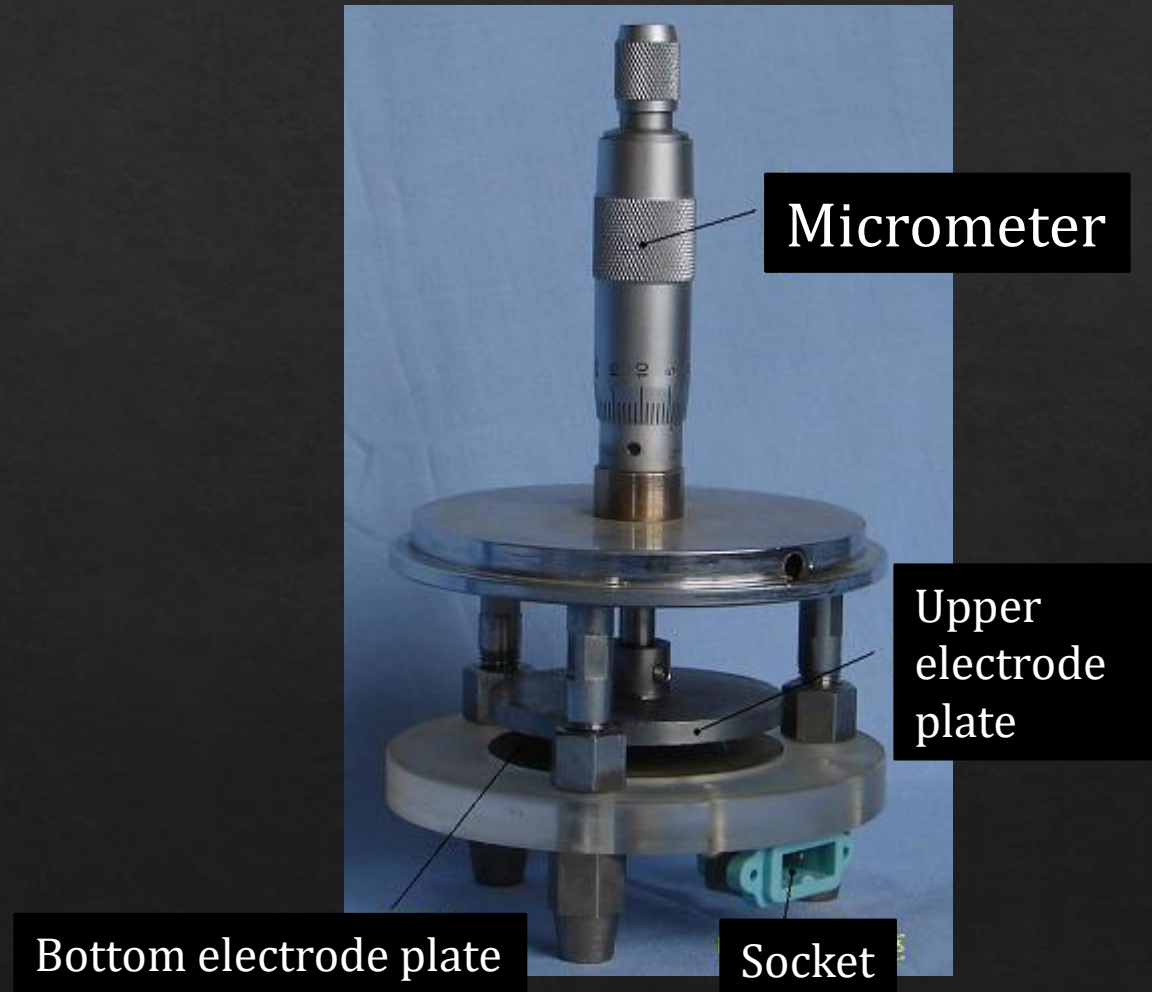

2、 Alternating current bridge

We choose 10kHz
and read from
'DISPLAY A' section.



3、 Solid dielectric sample

Solid samples are disk-shaped PTFE plastics, the material with high temperature resistance, corrosion resistance, high dielectric strength characteristics, therefore, commonly known as "king of plastic".



Steps and Contents:

1、 Measure the dielectric constant of solid dielectrics by alternating current bridge

- (a). Measure the diameters of the sample as d , three times using Vernier caliper. Measure the thickness of sample as t three times using micrometer.
- (b). Let the two plates touch each other, and read the starting value of the micrometer as D_0 .
- (c). Set the distance of the plates as 5.500mm. Note that the real value in micrometer should be $5.500\text{mm} + D_0$. Read the capacitance data from 'DISPLAY A' section as C_1 .
- (d). Insert the sample carefully, read the capacitance again as C_2 .
- (e). Repeat steps (d) and (e) twice.



2、 Measure the vacuum dielectric constant and capacitance

Remove the sample, set the plates distance to 1.000、
1.100 ... 1.900 mm, and
Read the corresponding capacitance C .

$$C = \frac{\epsilon_0 S_0}{D} + C_F$$

Set $x = 1/D$, $y = C$

Calculate C_F and ϵ_0 using origin software.

Data processing:

Table I

Starting value $D_0 = \underline{\hspace{1cm}}$ mm, Distance $D = \underline{5.500}$ mm。

$$\epsilon_0 = 8.85 \times 10^{-12} \text{F/m}$$

	D/mm	d/mm	t/mm	S/mm ²	C ₀ /pf	C ₁ /pf	C ₂ /pf	C _C /pf	ε _r
1	5.500								
2									
3									
Average									

d: diameter of sample,

t: thickness of sample,

$$S = \pi(d/2)^2$$

$$C_0 = \frac{\epsilon_0 S}{D}$$

$$C_C = C_2 - C_1 + C_0$$

$$\epsilon_r = \frac{C_C \cdot t}{\epsilon_0 S - C_C (D - t)}$$

Table II

D/mm	1.000	1.100	1.200	1.300	1.400	1.500	1.600	1.700	1.800	1.900
1/D(mm-1)										
C/pF										


$$C = \frac{\epsilon_0 S_0}{D} + C_F$$

Set $x = 1/D$, $y = C$,
calculate ϵ_0 using Origin
software.

The Slope is the value of
 $\epsilon_0 S_0$,

The Intercept is C_F .

Note: $S_0 = 21.61 \text{ cm}^2$.

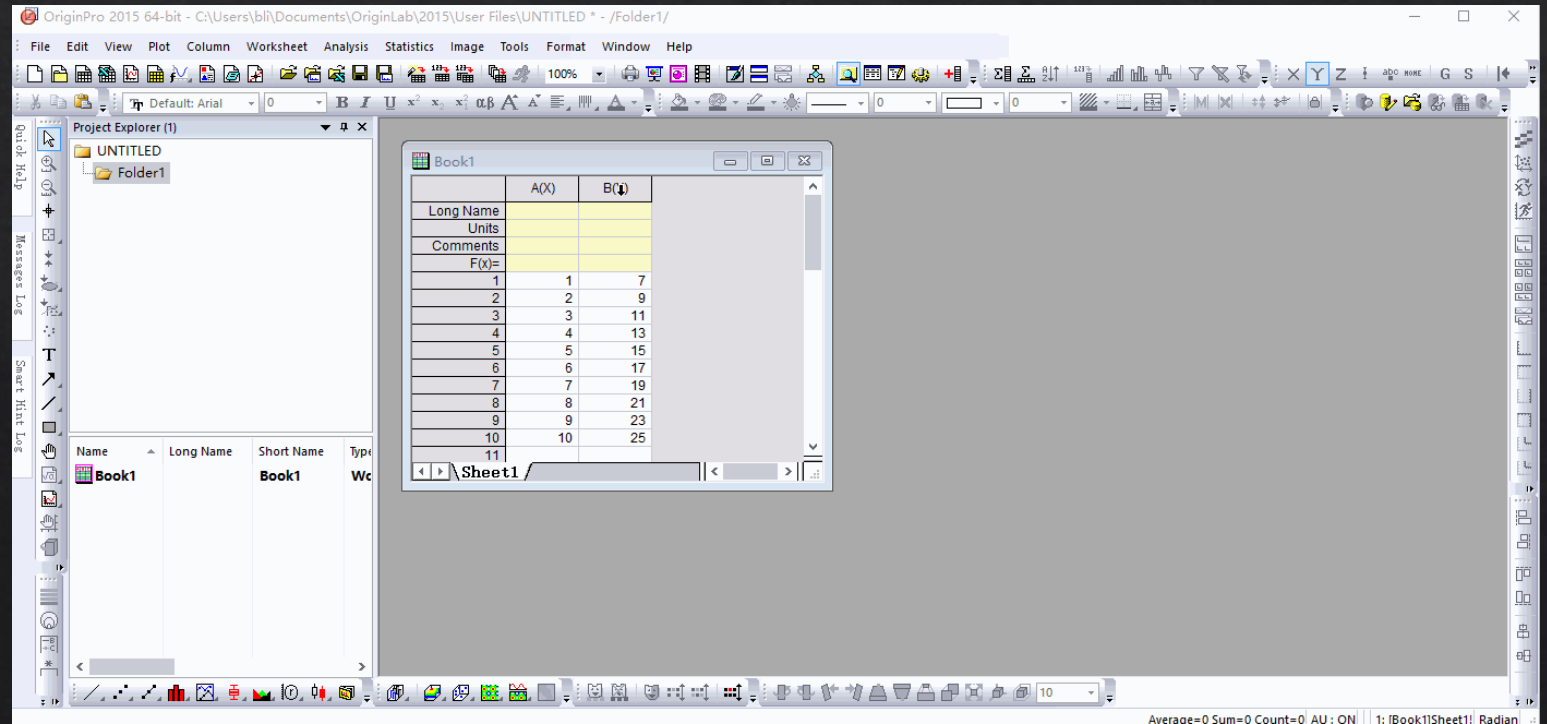

END